


## SL-5 EZ Flow Cell

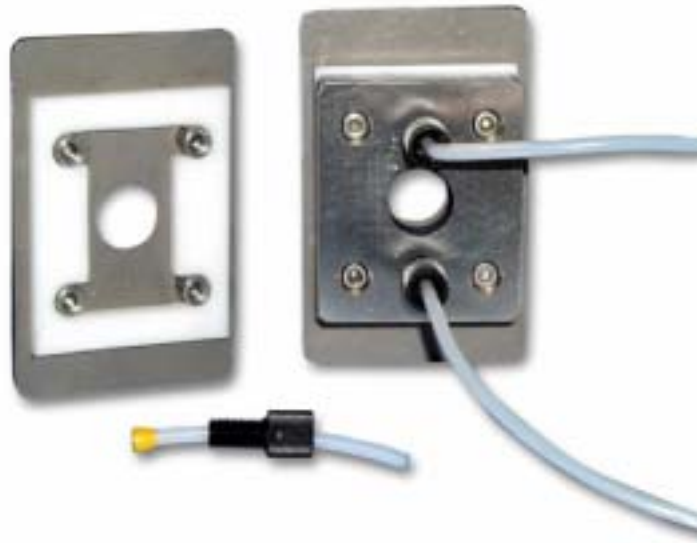

Our Model SL-5 E-Z™ Flow Cell is available in both sealed and demountable liquid flow cell configurations. The SL-5 E-Z™ Flow is made from stainless steel and features a Delrin® optical alignment cavity which, in the demountable version, precisely positions either a sandwich cell or the windows, spacers and gaskets in proper alignment with the cell filling ports. The cell features removable and replaceable compression nut and ferrule assemblies that hold the tubing in place. The standard nut is Delrin® and the standard ferrule is Tefzel®. However, stainless steel nut and ferrule assemblies are available. The cell can be used with both 1/16" OD and 1/8" OD tubing, and tubing diameters can be changed simply by removing the ferrule and nut assembly and inserting an inexpensive replacement. Special order high pressure versions of this cell are also available.

0007-7325 SL-5 Flow Cell, demountable with compression fittings for comprising a Delrin® compression nut sealed with a Tefzel® ferrule (specify 1/8" OD tubing or 1/16" OD tubing)

0007-XXXX Sealed SL-5 E-Z™ Flow Cell, same as 0007-7325 but as sealed cell

Please select 38.5\*19.5\*4 mm window for above